



Life CPAP Nasal Mask



- Soft dual-layer silicone cushion for better seal, comfort and stability to enhance patient compliance.
- Lightweight design, available in three sizes.
- Low flow resistance.
- 4.3mm port for CO₂ sampling, pressure monitoring or providing supplemental gas.

Features

For continuous improvement, product specifications subject to change without prior notice.

Lite CPAP Nasal Mask is designed to fit you comfortably.

1 Quiet

Vent port to reduce turbine and noise.



2 27° Ergonomic Angle

Ergonomic angle design directs gas away from the user and bed partner.

3 Multifunctional Port

Monitor CO₂ rebreathing, pressure or provide supplemental gases.

4 15mm Elbow

15mm swivel elbow with light weighted tubing.

5 Sliding Forehead Adjustment

Flexible adjustment of displacement and angle of forehead support.



6 Adjustable Forehead Support

Three types of forehead pads to choose from. Improve stability during sleep.



7 Ergo Seal

Air-tight dual layer silicone cushion. Reduce pressure sore.



8 Flexible

Easy to wear, quick release clip.



Order Information

REF	Descriptions	Packaging: 20 sets/box
ZN12	Lite CPAP Nasal Mask 15M with Economy Headgear, S	
ZN13	Lite CPAP Nasal Mask 15M with Economy Headgear, M	
ZN15	Lite CPAP Nasal Mask 15M with Economy Headgear, L	
ZN17	Lite CPAP Nasal Mask 15M with Economy Headgear, S	
ZN18	Lite CPAP Nasal Mask 15M with Economy Headgear, M	
ZN19	Lite CPAP Nasal Mask 15M with Economy Headgear, L	

DeZire®



Lite Nasal Mask

Made of lightweight, durable plastic material and adjustable forehead support, designed to fit comfortably and improve sleep quality. Ergo seal and forehead pad to reduce air leakage and eliminate pressure sore, maximizing therapy efficiency and compliance. Special vent port, designed to reduce air turbulence and ambience noise level, giving you and your partner a quieter night sleep.



Light and quiet,
3 sizes available



Vent port to reduce air
turbulence and noise



Quick release clip for
more flexibility



Multifunctional port for
monitoring or
supplemental gases



Soft dual-layer silicone
cushion for better seal,
comfort and stability



Slidable forehead
support pad

"40% lighter compared
to most common
models on the market."


[Product Feature]

1. 현존하는 Nasal Mask 중 가장 가벼운 무게의 제품입니다. 기존 제품 대비 40% 가볍습니다.
2. CPAP 운용 시 발생하는 소음을 최소화 하기 위한 'Vent Port'가 있어서 '조용한 마스크' 입니다.



3. Small/Medium/Large 의 3가지 사이즈가 있습니다. 보통 여성분이나 얼굴 크기가 상대적으로 작으신 분들이 Medium을 사용하고 남성 분의 경우 Large 사이즈를 선택하십니다.

4. Headgear(헤드밴드)와 마스크의 연결 및 탈착이 매우 간단합니다.

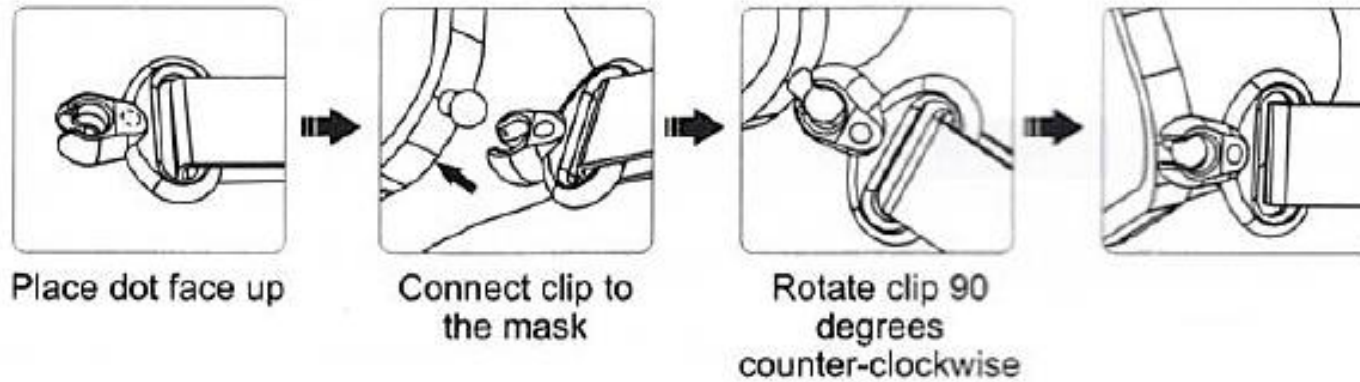

5. Multifunctional Port가 있어서 양압기(CPAP)의 부가적 기능인 CO2 sampling 및 Pressure Monitoring의 적용이 가능합니다.



6. Dual Layer의 매우 부드러운 인체친화적인 실리콘 재질의 쿠션이 있어 마스크 착용 시 코 주변부에 가하는 압박을 최소화하였으며 편안함을 증대시켰습니다.



Dual Layer의 매우 부드러운 인체친화적 실리콘 재질의 쿠션

7. 이마 지지대가 조절이 가능하여 마스크의 개인별 맞춤에 최적화된 마스크 착용이 가능합니다.

(여분의 이마 지지대 1개 포함.)



[제품 구성]

